



www.rhopointinstruments.com

TRANSMISSION HAZEMETER

RHOPOINT
NOVO-HAZE 

- Measures total transmission and haze
- Measure to ASTM D1003 (CIE C)
- Suitable for plastic films, and other transparent materials


Rhopoint Instruments Limited
is a member of the
Rhopoint Group

NEURTEK
i n s t r u m e n t s

T. 943 82 00 82 / info@neurtek.com
www.neurtek.com

Fast and accurate measurement of the optical quality of plastic films, and other transparent materials.

This instrument measures total transmission and haze according to ASTM D1003 (CIE C), the most important standard used in most QA applications.

Manufactured in direct response to industry requirements, the Novo-Haze TX is offered at huge savings compared to other instruments which contain additional superfluous test methods.

The instrument features an intuitive user friendly interface which minimises the test time and makes it an ideal choice for both QC and R&D.

Uncompromising design, high quality materials and European manufacture make the Novo-Haze TX the ideal choice for any laboratory or QA environment.

At 50% of the price of the market leader for only measuring to this standard (ASTM D1003, CIE C), the Rhopoint Novo-Haze TX represents huge savings without compromising measurement accuracy.



WHAT IS TRANSMISSION HAZE?

The Novo-Haze TX transmission hazemeter quantifies the optical qualities of transparent materials. The two most important aspects are:

1. Transmission – a measurement of the total amount of light that passes through a material that is influenced by absorption and reflective properties.

In the plastic film industry the transmission measurement is related to the opacity of the material and/or the hiding power of the applied coatings.


2. Transmission Haze – a measurement of the light scattering characteristics of a material. Haze can be due to suspended particles or contaminants within a sample or fine surface texture & contamination. Haze measurement can be used to quantify optical characteristics of plastics and packaging films.

In packaging applications hazy films may reduce a consumer's quality perception as packed products can look cloudy and indistinct. For plastics with haze the visibility of the test material becomes more apparent and the contrast of viewed objects is reduced.



Viewed through material with low haze


Viewed through material with higher haze

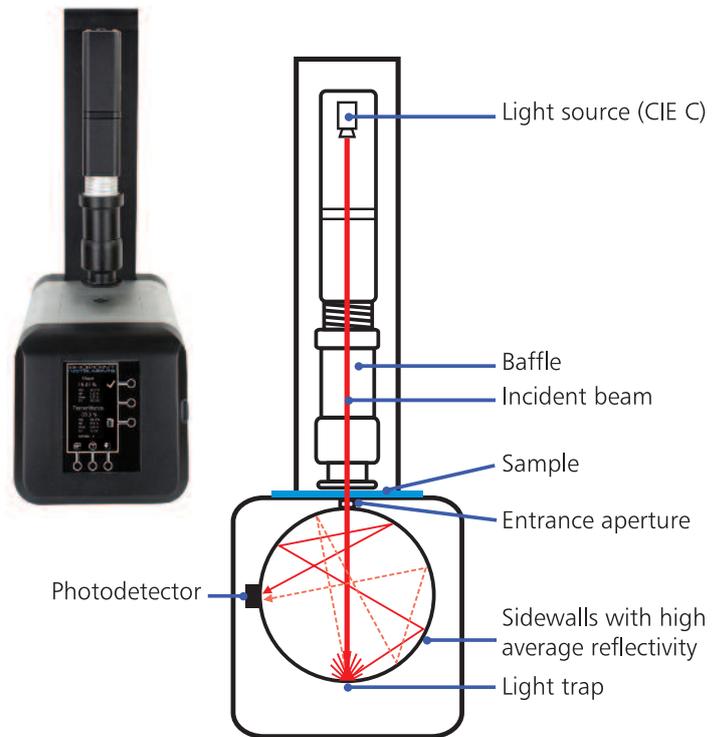

USER INTERFACE

The Novo-Haze TX has an intuitive touch sensitive interface making it incredibly easy to use. There are no complicated menu structures – all measurement functions are accessible from one screen.

Single touch button or the footswitch can be used to initiate single measurements or automatic measurement mode.

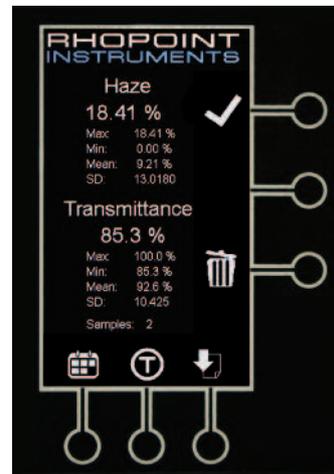

At the end of the test the percentage of Haze and Transmittance are shown on screen.

If multiple tests are performed, statistics for the batch will be displayed and can be printed to the results printer to keep with retained samples.



SPECIFICATIONS

HAZE

- **Repeatability**
± 0.05 units SD
- **Reproducibility**
± 0.4 units SD
- **Resolution**
0.01%

TRANSMISSION

- **Repeatability**
± 0.2 units SD
- **Reproducibility**
± 1.0 units SD
- **Resolution**
0.01%

Light source

- White light LED filtered for equivalence to CIE C luminosity function λ plus human photopic response

Measurement Area

- \varnothing 14mm

Measurement Range

- 0-100%

Geometry

- 0° / diffuse

Measurement Standard

- ASTM D1003 (CIE C)

Memory

- >10,000 via USB stick
- Internal memory: 99 (for statistical analysis only)

Power Supply

- 230V / 50Hz, 115V / 60Hz

Operating Conditions

- Operating temperature: 15 to 40°C
- Storage temperature: 0 to 50°C

Data Transfer

- Data transfer via USB memory stick to csv

Dimensions & Weight

- 580mm x 200mm x 320mm (H x W x D)
- Weight: 16.9kg
- Packed weight: 21.7kg
- Packed dimensions: 720mm x 440mm x 320mm (H x W x D)
- Commodity code: 9027 5000

Order Code

- Transmission Hazemeter: A3000-001

INCLUDED ACCESSORIES

- 2 x Haze checking standards, NPL traceable
- 24V DC power supply with mains cable
- Calibration certificate
- USB data stick

OPTIONAL ACCESSORIES

- Results printer



RHOPPOINT INSTRUMENTS www.rhopointinstruments.com (+44)1424 739622	
Transmission Hazemeter	
Instrument ID:	TN410018
Last Certified:	04/10/16
Test Date:	Oct 04 2016
Test Time:	Mon 09:40:42
Haze: 18.41 %	
Max: 18.41	Min: 0.00
Mean: 9.21	SD: 13.0180
Transmittance: 85.3 %	
Max: 100.0	Min: 85.3
Mean: 92.6	SD: 10.425
Samples: 2	

- Footswitch
- Custom sample holders

CALIBRATION AND SERVICE

Fast and economical service via our global network of accredited calibration and service centres. Please visit www.rhopointinstruments.com/support for detailed information.



Certificate no: FM 29741
ISO 9001:2008



LOCAL AGENT

NEURTEK

i n s t r u m e n t s

T. 943 82 00 82 / info@neurtek.com
www.neurtek.com

00564C/10/16